


# MYO/PUNCH ROTATING SURGICAL PUNCH

Designed to Cut a Clean and Round Hole

- Hollow cutting tip with 10 degree undercut moves tissue away from inside edge of hole during cut to provide a clean and round hole.
- Dual cutting action combined with a full 90 degree rotation of a scalpel sharp edge creates a clean circular hole.



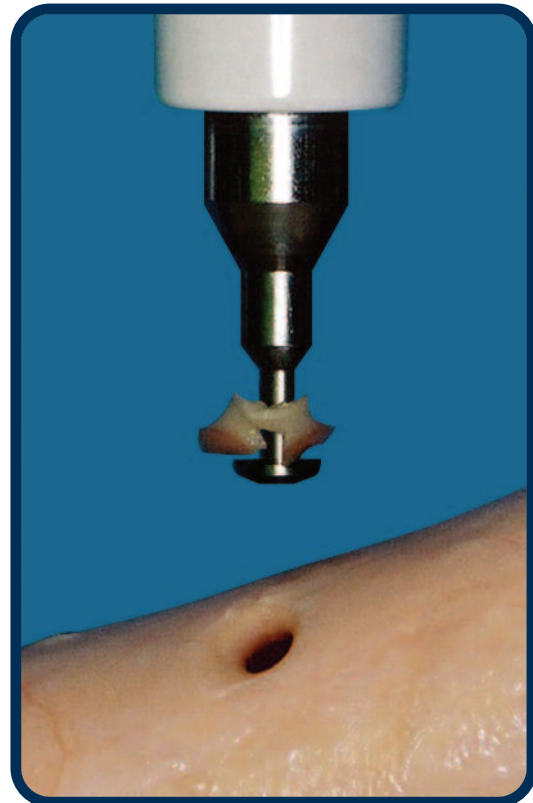


## Syringe Style Grip

Designed to provide comfort during use and excellent visibility of the operative site.

## Hollow Cutting Tip

Hollow cutting tip with 10 degree undercut moves tissue away from inside edge of hole during cut to provide a clean and round hole.

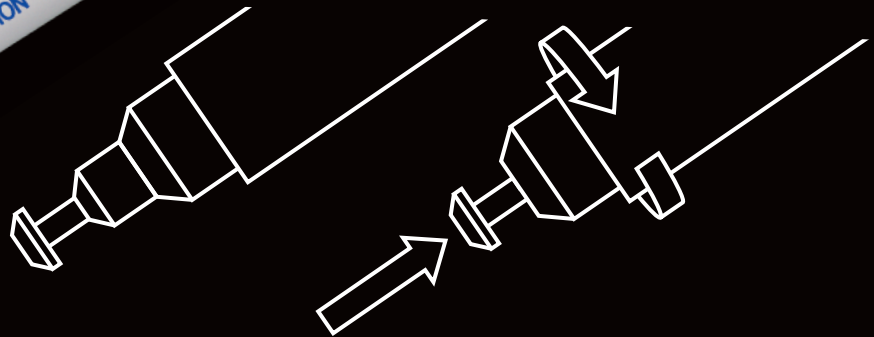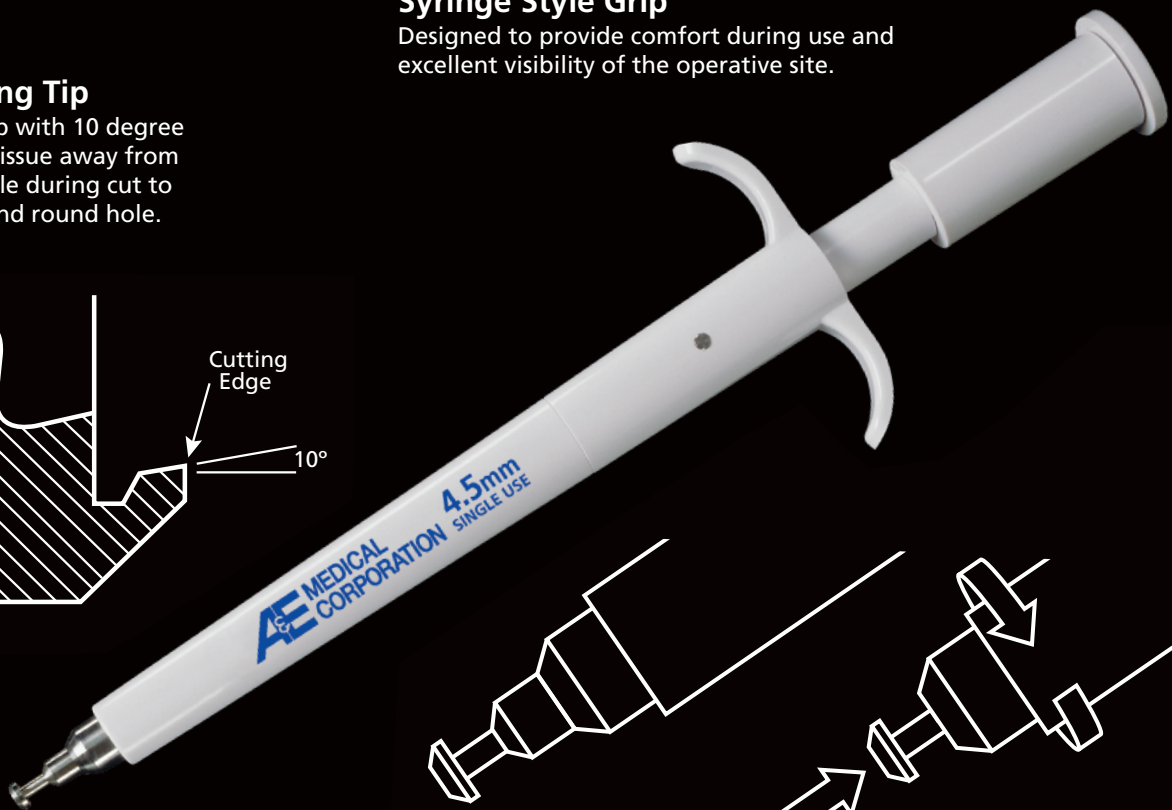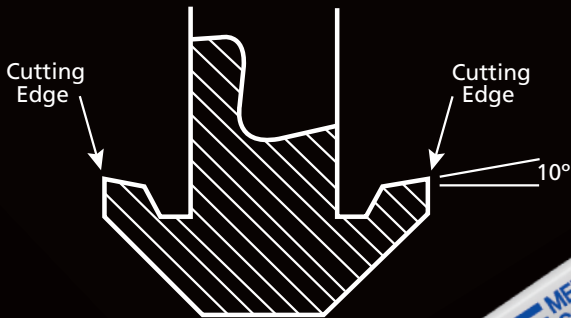

## Low Profile Cutting Tip

Designed for easier entry.

## Patented Rotating Cutting Action

Dual cutting action combined with a full 90 degree rotation of a scalpel sharp edge creates a clean circular hole.

## ORDERING INFORMATION

<u>CAT. No.</u>	<u>DESCRIPTION</u>	<u>DIAMETER</u>	<u>PACKAGING</u>
080-271	Rotating Surgical Punch	2.7mm	Sterile - 6 per box
080-351	Rotating Surgical Punch	3.5mm	Sterile - 6 per box
080-401	Rotating Surgical Punch	4.0mm	Sterile - 6 per box
080-451	Rotating Surgical Punch	4.5mm	Sterile - 6 per box
080-481	Rotating Surgical Punch	4.8mm	Sterile - 6 per box
080-501	Rotating Surgical Punch	5.0mm	Sterile - 6 per box

Patent Number 5,129,913